



ARJO HYGIENE AND WELLNESS SOLUTIONS

Therapeutic and tailored assisted bathing and showering

Assisted bathing and showering

Arjo integrated hygiene and wellness solutions enable you to deliver a tailored and therapeutic experience for residents and patients, helping you improve wellbeing and efficiency in your care environment.

Successfully meeting the hygiene needs of residents and patients in both acute and long-term care is not without significant challenges. Of all handling tasks, hygiene care is considered among the most high risk.

At Arjo, we believe bathing and showering should be deeply calming with sensory experiences that promote wellbeing and improve quality of life for residents and patients. Our solutions are specifically designed to help you deliver one-to-one care that is beneficial for your residents and patients, and safe and productive for caregivers.

- Carry out your daily hygiene tasks efficiently and ergonomically
- Promote a dignified experience that is calm and therapeutic
- Coverage across a wide spectrum of resident and patient needs
- Designed to be used by a single caregiver



The Mobility Gallery

Understanding the unique needs of individual residents and patients is key to quality hygiene and wellness care.

Our Mobility Gallery™ is an assessment tool that helps you find the right hygiene solution depending on a resident or patient's specific mobility levels as well as their personal preference, as indicated in the scale below.

Albert

Who is ambulatory and independent.



Barbara

Who can support herself to some degree.



Carl

Who sits in a wheelchair and has little capacity to support himself.





Doris

Who has no capacity to support herself.



Emma

Who is almost completely bedridden and totally dependent.



All Arjo hygiene and wellness solutions are categorized according to this scale, aligned to the resident and patient for which they are most suitable.

Assisted showering

Hygiene needs differ significantly across care settings, from patients that require basic assistance to residents with impaired cognitive and physical capabilities. Our assisted showering solutions allow you to work in an ergonomically sound position while maintaining eye contact with your resident or patient, supporting a beneficial and rewarding interaction.

No resident or patient is ever the same and needs and preferences may vary from day to day. Arjo showering solutions are designed to help you optimise the delivery of daily hygiene care by tailoring the showering experience to the unique needs of your resident or patient.

Carino

A comfortable, supportive shower chair for everyday assisted hygiene tasks such as washing and toileting, Carino® is height-adjustable to reduce static load and is also designed for intuitive ease of use. As the Carino is elevated, the resident or patient experiences a comfortable and secure reclined position



Designed with Barbara and Carl in mind






Carendo

A multipurpose hygiene chair, Carendo® helps reduce static load and contributes to a comfortable and dignified experience for your dependent residents or patients, even during the most sensitive hygiene tasks. The Carendo's unique Care Raiser™ function allows for safe and efficient undressing, dressing and incontinence care.

UNIQUE CARE RAISER FUNCTION



Designed with Carl and Doris in mind

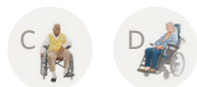

Carevo

An innovative solution for immobile residents and patients, the Carevo® shower trolley's patented design delivers ergonomic efficiency that helps you optimize the showering experience.

PRESSURE-RELIEVING FLEXI ZONE SURFACE

ERGO ACCESS AREA



Designed with Emma in mind



Assisted bathing and transfer solutions

Bathing promotes a sense of well-being for residents and patients and enhances quality of life, which is why our range of hygiene and wellness solutions are designed with flexibility in mind, ensuring all residents and patients can experience the therapeutic benefits of bathing. Our bathing solutions are fully integrated with the Arjo range of mobility and lifting aids, ensuring that you can safely and comfortably transport residents and patients from bed or wheelchair into the bath – and back again – in one seamless routine.

Recumbent bathing

To ensure a comfortable and safe bathing experience, Arjo has designed bath lifts that allow a single caregiver to transfer a resident or patient without the need of manual lifting.

You can also select your bath lift or other mobility device depending on the mobility level of your resident or patient and their preferred way of bathing (recumbent, reclining or sitting).

Our family of bath lifts is compatible with Arjo recumbent baths. Sling bathing via MaxiMove® and/or Arjo ceiling lift is also possible in Arjo's Parker® or recumbent bath range for more dependent residents' comfort. Sling bathing can be performed in any of our models, including Parker and Malibu™.

Alenti

Designed for active or semi-active residents and patients, the Alenti® bath lift chair provides a comfortable sitting position with the ability to recline as desired. The lift features a small footprint design, is easy to manoeuvre and can be docked from either side of the bath tub enabling a tailored bathing routine. Additional essential hygiene tasks, such as toileting, can be completed while your resident or patient is seated on the Alenti.



Designed with Albert, Barbara and Carl in mind





Miranti

Providing maximized powered support, the Miranti™ lift and bathing trolley ensures a single caregiver can position a dependent resident or patient into a reclining position for a beneficial bathing experience. Miranti's power drive function allows for effortless controlled transport from bed to bathing room and back.



Compatible with:
 Rhapsody® & Primo® 2,270 mm (89 3/8") length
 Rhapsody 2,520 mm (99 1/4") length
 Harmonie 2,050 mm (80 3/4") length
 Designed with Carl, Doris and Emma in mind



Maxi Sky 2

The Maxi Sky® 2 ceiling lift ensures transfers are fully ergonomic and free from strenuous handling. The ceiling lift with mesh sling allows water to pass through the material easily, making it an ideal choice for the showering and bathing of heavily dependent residents and patients.



Designed with Carl, Doris and Emma in mind



Recumbent bathing solutions

The Arjo recumbent bathing solutions enable an ergonomic and therapeutic bathing routine across the full spectrum of dependent residents and patients and are designed to enable easy maneuvering and docking of hygiene lifts where space is limited.

Arjo's recumbent bathing solutions include a range of unique bath shapes and is fully compatible with the range of Arjo hygiene lifters, giving you the flexibility to choose the right bathing solution for your resident or patient. All baths are fully ergonomic and height adjustable, ensuring you are always able to work in a safe and comfortable position and can enjoy a positive one-to-one interaction throughout your bathing task.

All Arjo recumbent baths include Hydromassage™, Hydrosound™, or Still Bath choices, along with optional Sound & Vision™ multi-sensory solution, designed to help you optimize the therapeutic benefits of the bathing experience for your resident or patient.

Rhapsody

A unique keyhole shape with additional shoulder space designed to improve access to your resident or patient. The solution provides greater access to the upper body and has been designed to enhance the comfort of the resident or patient.



Suitable for all mobility levels





Primo

The familiar shape of a traditional straight bathtub, Primo® provides efficient and comfortable bathing.



Harmonie

With a slight trapezoid shape, Harmonie is designed with a wider seat area to suit the needs of larger residents.



All bathing solutions include safety features:

- Thermoscopic mixing valve for pressure and temperature regulated water control
- Unique triple scalding protection with automatic shut off should water temperature exceed a safe level
- Integrated disinfectant dispensing system
- Emergency lowering

Suitable for all mobility levels



Height-adjustable reclining bathing solutions

With additional reclining functionality, the Arjo Parker is designed to provide more able-bodied residents and patients the therapeutic benefits of bathing.

Parker

The height-adjustable Parker reclining sit bath is designed to help you provide a safe and beneficial bathing experience that is easily accessible for your resident or patient through a specialized side-door entrance. Designed for residents in categories A to D, the door entry removes the need for an extra hygiene bath lift, and additions such as Air Spa and Sound & Vision offer a more pleasurable bathing solution.



Designed with Albert, Barbara, Carl and Doris in mind







Bathing and wellness

Therapeutic
benefits
that go beyond
daily routine

Bathing promotes a sense of wellbeing for residents and patients: a warm bath can help ease muscles, improve skin condition and relax the mind. The Arjo range of bathing solutions ensures residents and patients across the full spectrum of mobility can experience the therapeutic benefits of bathing.

Arjo bathing solutions incorporate a range of features and therapy options that enhance the bathing experience and assist you in creating an enjoyable as well as beneficial one-to-one interaction.




Hydromassage* — spa therapy

The built-in Hydromassage™ unit enables you to create whirling water with millions of small bubbles in the bath, providing a gentle massaging effect that can stimulate peripheral blood circulation, which may ease muscle tension, for a pleasant and sensory experience.



Hydrosound — cleaning therapy

Hydrosound helps to protect the skin of your residents and patients by offering an alternative to whirling water. Hydrosound sends sound pressure waves through the water, creating millions of microscopic bubbles that implode and act as gentle "skin scrubbers," removing dirt particles from the skin without invasive manual washing.

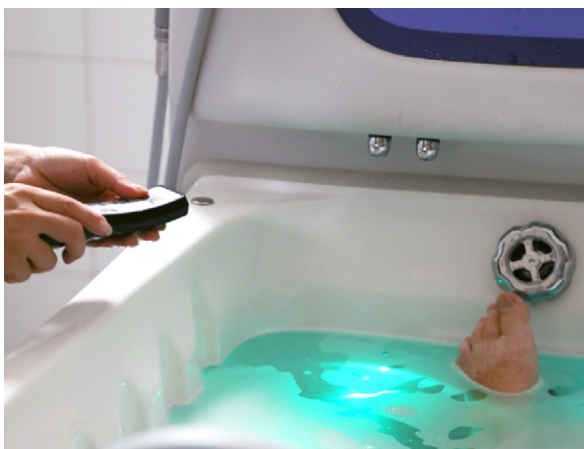

Sound & Vision** — multisensory experience

Sound & Vision provides pleasant, colorful underwater lighting and soft music that enhance the power of bathing to help calm stressed or anxious residents and patients. Underwater lights, mounted at the foot of the bath, can be set to rotate between red, green, blue, pink and yellow, while a USB stick with MP3 files enables you to play soft and relaxing MusiCure® tracks or music to the liking and benefit of both dependent and cognitively impaired residents and patients.



Air Spa

Providing residents and patients with a gentle and relaxing spa sensation, Air Spa can be adjusted from a low level setting to a more invigorating massage sensation.

* Only for Arjo recumbent baths

** Only Arjo Parker and Arjo recumbent baths

Overview



Carino		●	●		
Carendo			●	●	
Carevo					●
Alenti	●	●	●		
Miranti			●	●	●
Maxi Sky			●	●	●
System 2000 Rhapsody	●	●	●	●	●
System 2000 Primo	●	●	●	●	●
System 2000 Harmonie	●	●	●	●	●
Parker	●	●	●	●	

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure injuries and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind. www.arjo.com

Arjo Inc. · 2349 West Lake Street · Addison, Illinois 60101 · 800-323-1245 · www.arjo.us

Arjo Canada Inc. · 90 Matheson Blvd West, Suite 350 · Mississauga, ON, L5R 3R3 · 800-665-4831 · Info.Canada@arjo.com · www.arjo.ca