

Gravity 8 Deluxe Long Term Care Pressure Redistribution Mattress

HCPCS:

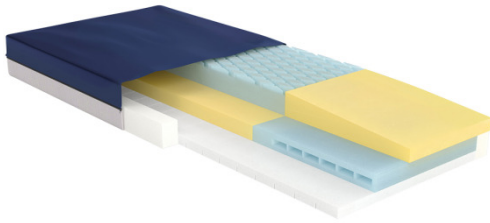


FIGURE A

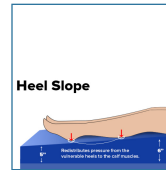


FIGURE B

- High quality therapeutic foam mattress provides pressure redistribution and shear/friction reduction
- High end therapeutic foam mattress maximizes value and comfort with 3 full layers of latex free foam
- Full width top layer uses unique die-cutting to optimize zoned pressure redistribution
- Bottom layer provides full length horizontally scored articulation cuts extending the durability and life of the mattress
- High density foam perimeters provide solid edge support and facilitate safe patient transfers
- Available with or without raised side rails
- Heel section is sloped with channel cut foam center offloading pressure from vulnerable heels
- Removable multi-stretch cover is fluid resistant, low shear and vapor permeable helping protect the skin from friction and moisture
- Meets flammability standard CFR 1632 and CFR 1633

PRODUCTS	Description
15870	Gravity 8 Mattress 80", 1/EA
15884	gravity 8 Mattress 84", 1/EA
15886	Gravity 8 Mattress 76", 1/EA
15888	Gravity 8 MattressW/Elev Perim 80", 1/EA

SPECIFICATIONS	Description
Warranty	10 Year Limited Foam 3 Year Cover
Weight Capacity	400 lbs
FireRating	CFR 16 Part 1633
Side Rails	Yes
Cover Material	Stretch, Fluid-Resistant Nylon Top Waterproof Vinyl Bottom
Mattress Layers of Foam	3

ACCESSORIES	Description
15870C	Mattress Cover 15870C,1/ea, EA
15884C	Mattress Cover 15884C,1/ea, EA
15888C	Mattress Cover 15888C,1/ea, EA
15886C	Mattress Cover 15886C,1/ea, EA