


Sömn & Koppla

Customizable and layered foam mattresses

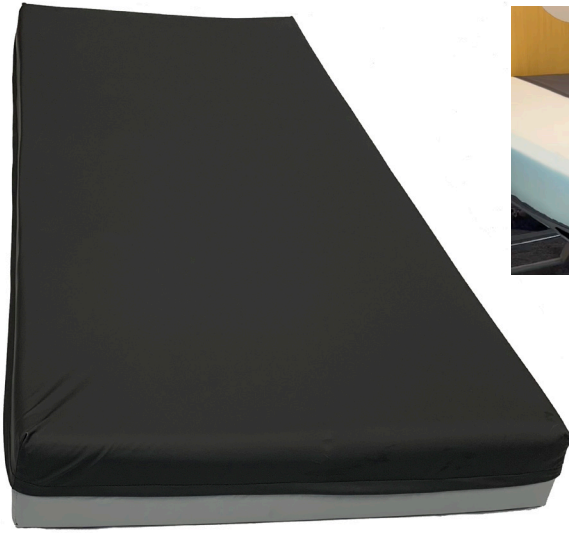

At a glance

- Ideal for long-term care clients
- 3" wide transfer edges to help reduce the risk of falls
- Fluid-resistant Startex full enclosure cover
- Sömn is firmer with a 4" base layer and 2" top layer with Serene™ foam at the heel
- Koppla is softer with a 2" thick foam base, plus a 2" middle, and a 2" top layer of foam with decreasing density
- Maximum client weight of up to 350-500 lb
- 4 widths: 36" / 39" / 42" / 48" x 80"L x 6"H
- 24-month warranty

The **Sömn** mattress features a dual-layered high resiliency foam for pressure redistribution and is the ideal choice for those looking for a firm mattress. It has a 4" thick base layer and a 2" top layer with lower density Serene foam at the heel for better immersion. Three-inch wide transfer edges extend to the top on both sides of the surface (shown above) to help reduce the risk of falls for long-term care clients.

The **Koppla** mattress features a multi-layered high resiliency foam for improved support and immersion. It has a 2" thick base layer, 2" middle, and a 2" top layer with decreasing foam density. Subtle transfer edges that are 3" wide go up the first 4" on both sides of the surface for added support. Soft Serene foam is used for the entire 2" top layer providing comfort and immersion and makes this the ideal choice for lightweight users and children.