


# Symmetry

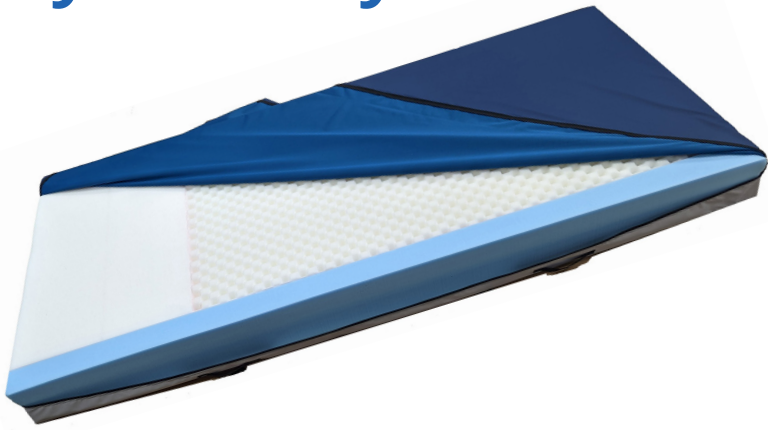

## At a glance

- Ideal for long-term care clients
- Wide transfer edges to help reduce the risk of falls
- Fluid-resistant, multi-way stretch cover
- Maximum client weight of up to 350-500 lb
- 4 widths: 36" / 39" / 42" / 48" x 80"L x 6"H
- 24-month warranty

Symmetry has three layers of foam with a convoluted high resiliency foam on top to improve client mobility, including a softer foam at the heel. It also includes firm side borders for secure transfers in/out of bed, and for keeping the client in a central position. Symmetry has a multi-way stretch, fluid-resistant cover which provides incontinent protection.